# aquapick®

CORDLESS ORAL IRRIGATOR

**aquapick®**  
CORDLESS ORAL IRRIGATOR

## Aquapick Co., Ltd

Адрес: 20F, SK V1 Center, 171, Gasan digital 1-ro, Geumcheon –gu,  
Сеул, Южная Корея  
Тел. 1588-2809 / E-mail. info@aquapick.ru

## User Guide

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# IMPORTANT SAFE GUARDS

All information in this section is provided to instruct you in the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

READ ALL INSTRUCTIONS BEFORE USING.

## 1. DANGERS

- ① Toddlers, infants, and people who are unable to operate should not be allowed to use this appliance.
- ② Never alter the appliance. Also, do not try to disassemble or repair it.
- ③ Do not place or store the product where it can fall or be pulled into tub or sink.
- ④ Do not use the product while bathing.

## 2. WARNINGS

- ① Never use the adaptor if it has a damaged cord or plug.
- ② Keep the cord away from heated surfaces.
- ③ The supply cord cannot be replaced. If the cord is damaged the appliance should be scrapped.
- ④ Please do not use it if you have lost the jet tip.
- ⑤ Do not fill the tank with water warmer than 40°C.
- ⑥ Do not submerge the device in water.

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- ⑦ Do not use for any purpose other than oral cleaning.
- ⑧ Keep the jet tip away from small children. The jet tip is small enough to be swallowed and can be lodged in the throat.
- ⑨ People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- ⑩ Do not press the jet tip too firmly against your teeth or gums. Doing so may injure your gums.
- ⑪ This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities or lack of experience and knowledge unless they have been given adequate supervision or instruction concerning use of the product by an adult responsible for their safety.
- ⑫ The appliance is only to be used with the power supply unit provided with the appliance.
- ⑬ Do not attempt to remove the battery.
- ⑭ The battery must be removed from the appliance before it is scrapped.

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
### MAIN UNIT

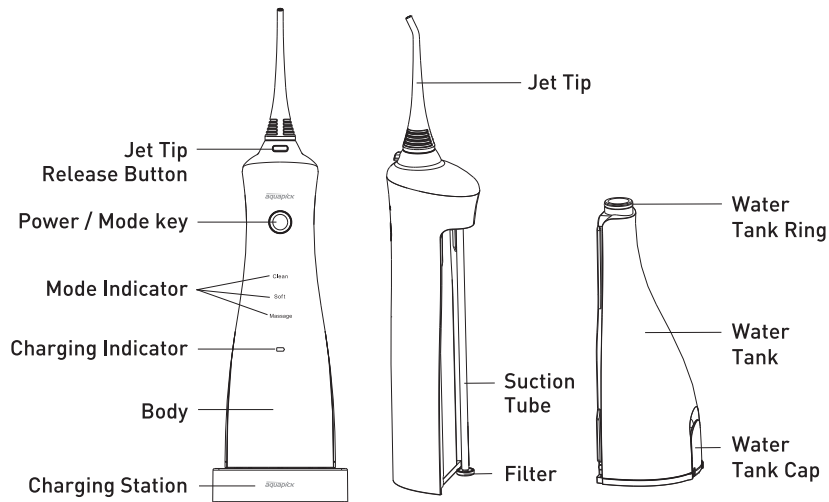
### JET TIP (2pcs)

Purchasable Separately  
1SET(2pcs/set)

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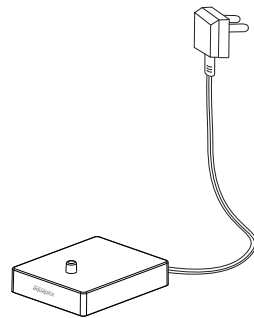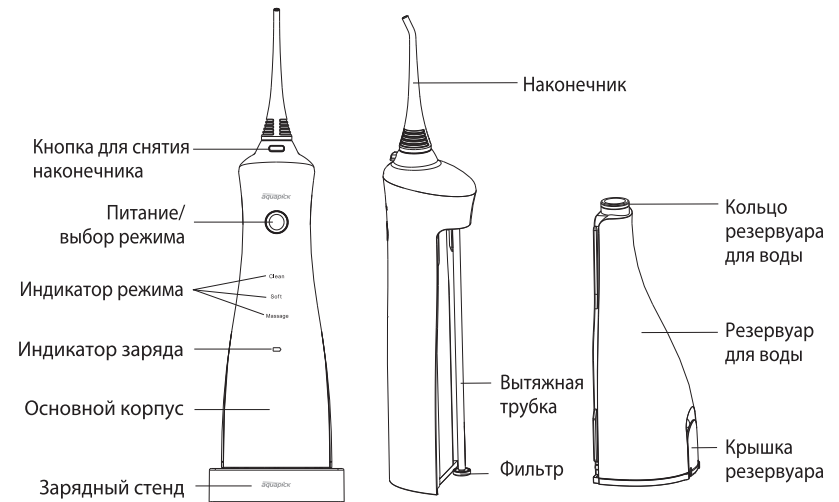
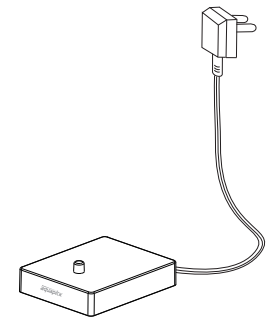
## GETTING STARTED

### Notes:

There might be some liquid remaining inside the main appliance before use (tank, water suction tube, etc), this is the remaining distilled water that was used for the performance (water flow) test, so there are no health concerns.


### Charging the device (Fig. 1)

- 1 Insert the adaptor plug into the power outlet.
- 2 Put the device up onto right place of the charger base.
- 3 The charge indicator blinks when the device is charging, and it will be steadily on after the device is full charged .

### Notes:

If the charge indicator is blinking during or after an irrigating, that means the battery need to be recharged.

It takes about 24 hours to fully charge it.

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


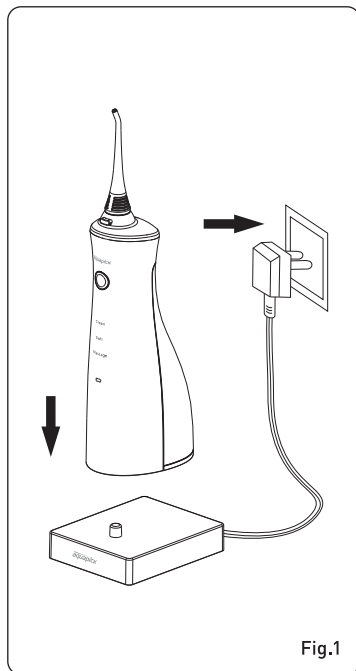
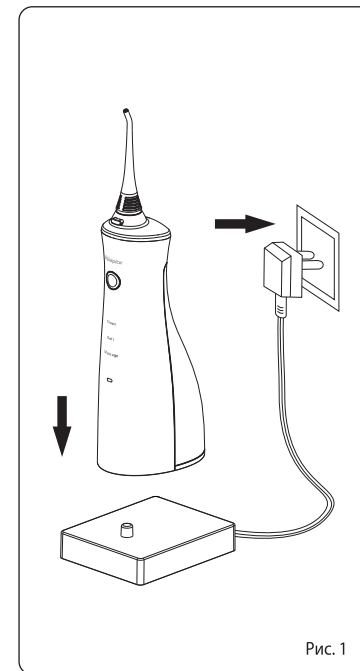
## HOW TO USE

### Use the oral irrigator

- 1 Attach the jet tip onto the unit. (Fig. 1)
- 2 Open the water tank cap to fill the tank with water whose temperature does not exceed 40°C, and then close the cap. (Fig. 2)
- 3 Or slide the water tank down off the unit, fill the tank with water whose temperature does not exceed 40°C, and then attach the tank onto the unit. (Fig. 3)
- 4 Press the Power/Mode key to turn on the oral irrigator. It starts working at Clean mode. During irrigating, if Power/ Mode key is pressed, it will work next mode and the mode indicator will light up. (Fig. 4)
  - a. Clean: for normal cleaning.
  - b. Soft: for gentle cleaning of sensitive gums.
  - c. Massage: for cleaning and massaging the gums.

After pressing the Power/Mode button at first, if the button is pressed within 5 seconds again, its mode is changed and if the button is pressed after 5 seconds, the power is turned off.
- 5 Hold the unit in an upright position with the jet tip point toward your teeth.
- 6 Keep your mouth opened slightly so that the water can flow out.

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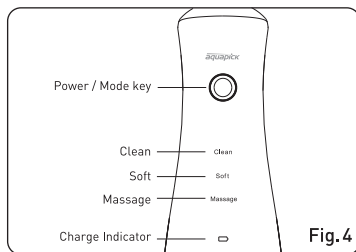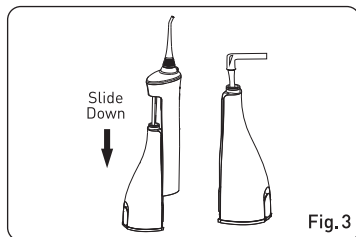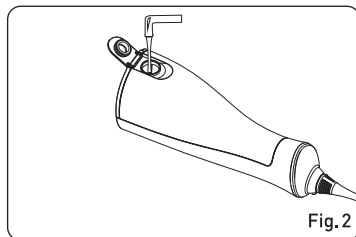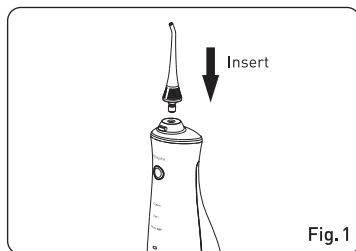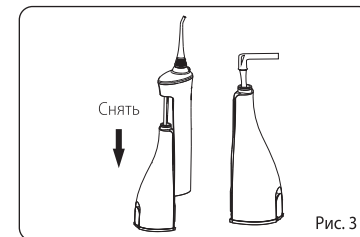
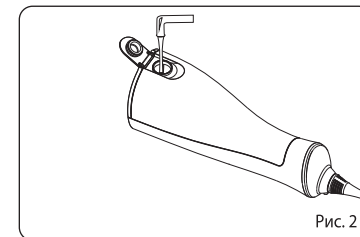
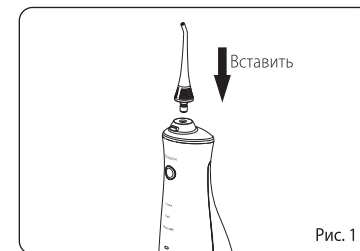
## HOW TO USE

- ⑦ Press the Power/Mode key to start irrigating. (Fig. 5)

- 1) Directing the water jet flow at a 90-degree angle to the teeth and gums.
- 2) Direct the water jet slowly along the teeth.
- 3) Clean between the teeth and around any braces, crowns or bridgework thoroughly.

- ⑧ Press Power/Mode key and then press the jet tip release button to replace other jet tip when you need. (Fig. 6)

- ⑨ When you are finished using this product, press Power/Mode key again to turn it off.

### After using the appliance (Fig. 7)

- ① Open the water tank cap and dispose of the water in the tank.
- ② Press Power/Mode key to turn on it, and drain the water inside the product.
- ③ Press Power/Mode key again to turn off the product, and close the water tank cap.
- ④ Wipe the product with a dry cloth.

### The remaining charge indicator

#### When charging




- ① Below 10%: Red lamp flickers
- ② Over 10%: Green lamp flickers.
- ③ Full charge: Green lamp always on.

#### When operating

- ① Below 10%: Red lamp always on.
- ② Over 10%: Green lamp always on.

#### Notes:

If you are not going to use this product for a long time, be sure to wipe and dry the product before storing it.

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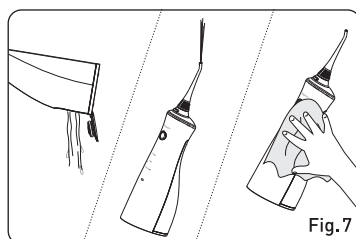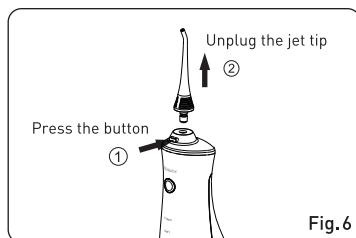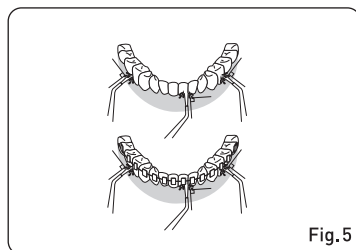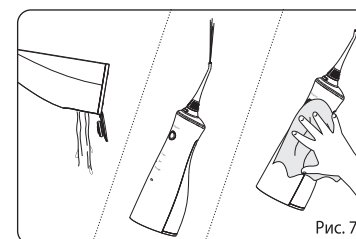
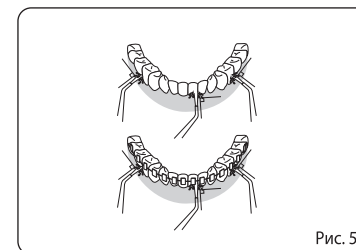
## MAINTENANCE

Use only water or neutral cleaning agent to clean the device. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may damage the device. Do not wash using hot water of over 50°C. Keep the device away from high temperature area and avoid direct sunlight.

### Main Unit

- 1 Remove the dirt by wiping using a wet cloth and wipe off using a dry cloth.
- 2 Do not submerge the main unit in water to clean.
- 3 Dry the product after using it.

### Induction Charger

- 1 Check the charger plug and cable at regular intervals, and wipe the adapter plug with a cloth.

### Jet Tip

- 1 Wash with water and wipe with a soft cloth.
- 2 We recommend replacing jet tips every 6 months for hygienic use.

### Water Suction Tube

- 1 Wash with water and wipe with a soft cloth.
- 2 Do not bend, pull or twist the water suction tube.

### Water Tank

- 1 Wash with water and wipe with a soft cloth.
- 2 Dry the tank when not using for a long period.

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




## DISPOSAL

The rechargeable battery inside your this product is easy to remove for recycling.

At the end of the product's life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below.

During below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

1. Remove two rubber feet on the bottom of the handle. (fig.8)
  2. Insert a screwdriver into the slot located at the bottom of the handle and turn counter-clockwise to release the bottom cap. (fig.8)
  3. Slide the water down off the device and unscrew the two screws beside of the suction tube. (fig.9)
  4. Open the top cover and put the inner out of the unit, then unscrew all the screws of the inner. (fig.10)
  5. Cut off the battery line and pull out the batteries from the inner. (fig.11)
- If you have questions about recycling, please contact your local waste management office.

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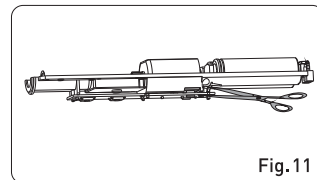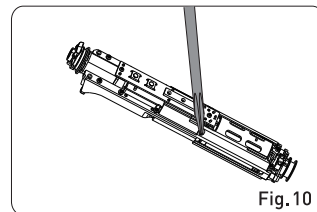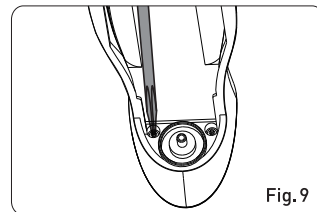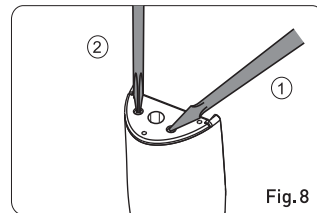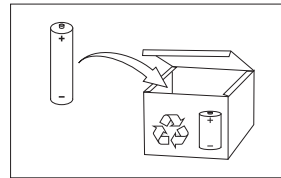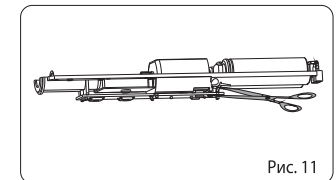
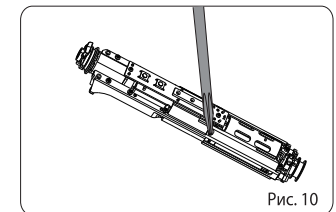
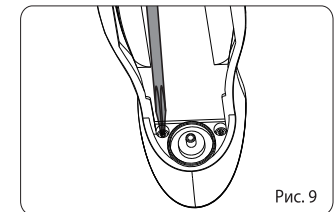
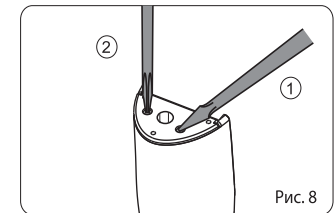
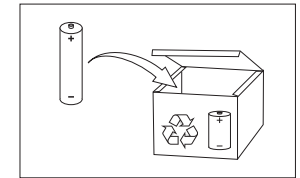
## TROUBLE SHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
It does not work	You've just purchased the product or it has been left unused for more than 3 months	Charge the unit for about 24 hours
It only works for a few minutes even after charging	The battery has reached the end of its lift cycle	Contact service center for support
	Charging time was insufficient	Charge the unit for about 24 hours
Water pressure is very low	There is dirt stuck in the jet tip	Replace the jet tip
	The jet tip is deformed	
	There is dirt stuck in the filter	Clean the filter or purchasable separately
Water does not come out	The tank is empty	Fill up tank
	The product is tilted too much while in use	Use with the product in a vertical state

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This product has passed strict testing and inspection. Aquapick Co., Ltd. guarantees quality for 12 months from date of purchase.

We will repair defective products that occur under normal use within the warranty period free of charge, and replace the product with a new one if it is difficult to repair.

### Precautions

Even during the warranty period, the following cases will be repaired at your expense:

- 1) The product is damaged due to improper use, mishandling or negligence after purchase.
- 2) The product is altered, modified or repaired by a party not authorized by Aquapick Co., Ltd.
- 3) The product is damaged as a result of connection to irregular voltage sources.

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